

ARTER & HADDEN, LLP  
1801 K Street, N. W. Suite 400K  
Washington, D. C. 20006-1301

JC553 U.S. PTO  
09/373984  
08/16/99

TO: Assistant Commissioner for Patents  
Box Patent Applications  
Washington D.C. 20231

Attorney Docket No. 70862/93137

UTILITY PATENT APPLICATION TRANSMITTAL

Transmitted herewith for filing is the patent application of:

**INVENTOR(S): Xing SU and Rolfe C. ANDERSON**

**TITLE: SINGLE-PHASE AMPLIFICATION OF NUCLEIC ACIDS**

In connection with this application, the following are enclosed:

APPLICATION ELEMENTS:

- ☒ Specification 25 pages; Abstract 1 page; Number of Claims 24  
☒ Drawings Total Sheets 1  
\_\_\_\_ Unexecuted Declaration and Power of Attorney - Total Sheets  
☒ Newly executed (original) Declaration and Power of Attorney  
\_\_\_\_ Microfiche Computer Program (Appendix)  
\_\_\_\_ Nucleotide and/or Amino Acid Sequence Submission (if applicable, all necessary)  
\_\_\_\_ Computer Readable Copy  
\_\_\_\_ Paper Copy (identical to computer copy)  
\_\_\_\_ Statement verifying identify of above copies

ACCOMPANYING APPLICATION PARTS

- ☒ Assignment Papers (cover sheet & document(s))  
\_\_\_\_ 37 CFR 3.73(b) Statement (when there is an assignee)  
\_\_\_\_ English Translation Document (if applicable)  
\_\_\_\_ Information Disclosure Statement (IDS) with PTO-1449. \_\_\_\_ Copies of IDS Citations  
\_\_\_\_ Preliminary Amendment  
☒ Return Receipt Postcard (MPEP 503)  
\_\_\_\_ Small Entity Statement(s)  
\_\_\_\_ Statement filed in prior application, status still proper and desired.  
\_\_\_\_ Certified Copy of Priority Document(s) with Claim of Priority  
(if foreign priority is claimed).  
\_\_\_\_ OTHER:

If a **CONTINUING APPLICATION**, check appropriate box and supply the requisite information:

☐ Continuation ☐ Divisional ☐ Continuation-in-part (CIP)  
of prior application Serial No.

☐ Amend the specification by inserting before the first line the following sentence: --This application is a ☐ continuation, ☐ divisional or ☐ continuation-in-part of application Serial No. , filed --

**CORRESPONDENCE ADDRESS:**

Don J. Peltó, Esquire  
1801 K Street, N.W. Suite 400K  
Washington, D. C. 20006-1301  
Phone: (202)-775-7980  
Fax: (202) 857-0172

**FEE CALCULATIONS:** (Small entity fees indicated in parentheses.)

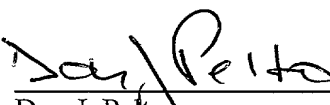
(1) For	(2) Number Filed	(3) Number Extra	(4) Rate	(5) Basic Fee \$760 (\$380)
Total Claims	24-20 =	4	x\$18 (x \$9)	\$ 72.00
Independent Claims	1-3 =	0	x \$78 (x \$39)	-0-
Multiple Dependent Claims			\$260 (\$130)	-0-
Assignment Recording Fee per property			\$40	\$ 40.00
Surcharge Under 37 C.F.R. 1.16(e)			\$130 (\$65)	-0-
			TOTAL FEE:	\$ 872.00

**METHOD OF PAYMENT:**

A check in the amount of \$ 872.00 is enclosed. This amount is believed to be correct. However, the Commissioner is hereby authorized to charge any deficiency or credit any overpayment to Deposit Account No. 50-0767.

Respectfully submitted,

Date: 16 August, 1999  
Docket No. 70862/93137

  
\_\_\_\_\_  
Don J. Peltó  
Reg. No. 33,754

5                   **SINGLE-PHASE AMPLIFICATION OF NUCLEIC ACIDS**

**CROSS REFERENCE TO RELATED APPLICATIONS**

This application claims priority to U.S. provisional application serial number 60/146,142, filed 30 July 1999.

10

**FIELD OF THE INVENTION**

The present invention relates generally to the amplification of nucleic acids in a single phase. In a preferred embodiment, the present invention can facilitate the proportional amplification of minute sample amounts of nucleic acids in a manner that may preserve the relative abundance of the individual nucleic acid species, or portions thereof, existing in the original sample.

15

**BACKGROUND OF THE INVENTION**

The isolation, characterization and manipulation of nucleic acids has numerous present or potential applications, including those in the basic research, diagnostic and forensic fields. Valuable information about gene expression in *in vivo*, *in situ*, and *in vitro* systems can be obtained by monitoring the abundance of the mRNA encoded by those genes. Methods involving the synthesis of cDNA from mRNA have also enhanced the study of gene expression, for example, by facilitating gene cloning and the production of desired recombinant proteins.

20

25

With existing methods for the study or use of mRNA and cDNA, one problematic scenario can arise where the sample size is small, or the relative abundance of an individual mRNA or cDNA species in a sample is low. In such situations, where the availability or accessibility of the desired mRNA or cDNA is compromised (or their amounts are otherwise limited), the lower limits of monitoring or manipulation systems may be exceeded, thus leaving the desired mRNA or cDNA undetected, unrecoverable or unworkable. Therefore, the amplification of such mRNA and cDNA is an important molecular biology methodology, with particular

30

0937394-0169

significance in facilitating the detection and study of a broader range of mRNA molecules, and the isolation and manipulation of mRNA available in only minute quantities.

Known nucleic acid amplification methods may typically involve multiple steps and varying reaction conditions such as organic extraction and precipitation. As result, these methodologies can be labor intensive and costly. For example, to amplify RNA in a heterogenous RNA population, known methods can require that a second strand cDNA synthesis is performed in a volume that is more than three times the volume of the first strand cDNA synthesis reaction. Prior to the *in vitro* transcription (IVT) reaction, in which RNA polymerase may use promoter containing double-stranded cDNA as templates for RNA transcription, the cDNA sample typically must be cleaned by, for example, phenol extraction and ethanol precipitation. When these cleaning and concentration steps are not used, the IVT reaction may be inhibited due to undesirable buffer conditions and enzymatic activities.

In addition, although methods exist for the amplification of nucleic acids, they generally suffer from a phenomenon known as biased amplification. In these cases, the amplified population does not proportionally represent the population of nucleic acid species existing in the original sample. This drawback may preclude meaningful or reliable conclusions regarding the absolute amount or relative abundance of a desired nucleic acid species in the tested sample.

One common problem encountered by past amplification methods is the preference for the amplification of shorter nucleic acid templates. The enzymes responsible for the production of complements or copies of the nucleic acid templates (e.g., DNA and RNA polymerases, or reverse transcriptases) achieve such synthesis through a sequential, oriented process, whether 5' to 3' or 3' to 5'. The probability that such an enzyme will complete a copying event thus may be greater with nucleic acid templates of shorter length. Accordingly, in a sample population containing nucleic acid templates of variable lengths, longer templates may be less likely than shorter templates to be amplified in complete, full-length form. This can result in a bias in the amplified population in favor of nucleotide sequences proximal to the 3' poly(A) tail of mRNA, for example, a phenomenon known as 3'-sequence bias.

5 The synthesis of longer templates can also be difficult or less efficient due to interference from secondary and tertiary structure in the template. For example, with respect to nucleic acid amplification based on polymerase chain reaction (PCR) methodologies, longer templates in a sample may be under-represented in the amplified product if respective primers cannot anneal to begin another round of copying because the first round did not proceed to completion. Other potential sources of bias can reflect relative differences between longer and shorter templates. For example, longer templates may (i) not denature sufficiently, or (ii) have a greater likelihood of mismatches, and thus error propagation through amplification, but (iii) have an ability to anneal more easily.

The foregoing shows a need for methods and products involving the amplification of nucleic acids in a simplified manner, and preferably, to facilitate the preservation of the relative abundance of the individual nucleic acid species existing in the original sample.

#### SUMMARY OF THE INVENTION

An objective of the present invention is therefore the amplification, preferably proportionally, of nucleic acids in a single phase.

15 In accomplishing these and other objectives, the present invention provides methods for the amplification, preferably proportionally, of nucleic acids that may comprise synthesizing double-stranded DNA from a single-stranded DNA population, and producing multiple copies of RNA from the double-stranded DNA, where the amplification occurs in a single phase. In one preferred embodiment, this process may incorporate the use of an automated machine, preferably a PCR thermocycler.

20 In a preferred embodiment, the present invention provides methods in which the single-stranded or double-stranded DNA population may be produced from a nucleic acid population selected from the group consisting of one or more of the following: genomic DNA, cDNA, total RNA, poly(A)<sup>+</sup> RNA, and oligonucleotides. In a preferred embodiment, the poly(A)<sup>+</sup> RNA may be mRNA.

25 The present invention also preferably provides methods, which may further comprise contacting the multiple copies of RNA with a solid support comprising nucleic acid probes, and detecting the presence or absence of hybridization of the

RNA to the nucleic acid probes on the solid support. In a preferred embodiment, the solid support, which may comprise nucleic acid probes, can be selected from the group consisting of a nucleic acid probe array, a membrane blot, a microwell, a bead, and a sample tube.

5           The present invention may preferably provide methods wherein the nucleic acid may be isolated from an eukaryotic cell or tissue, mammalian cell or tissue, or human cell or tissue. In a preferred embodiment, the nucleic acid may be isolated from a source selected from the group consisting of dissected tissue, microdissected tissue, a tissue subregion, a tissue biopsy sample, a cell sorted population, a cell  
10 culture, and a single cell. In another preferred embodiment, the nucleic acid may be isolated from a cell or tissue source selected from the group consisting of brain, liver, heart, kidney, lung, spleen, retina, bone, lymph node, endocrine gland, reproductive organ, blood, nerve, vascular tissue, and olfactory epithelium. In yet another preferred embodiment, the nucleic acid may be isolated from a cell or tissue source  
15 selected from the group consisting of embryonic and tumorigenic.

          In a preferred embodiment, the present invention may provide an amplified nucleic acid preparation comprising RNA obtained by the described methods. In another preferred embodiment, the present invention may also provide an amplified nucleic acid preparation comprising DNA obtained by the described methods.

20           The present invention preferably provides a gene expression monitoring system comprising a solid support, which comprises nucleic acid probes and the amplified nucleic acid preparations. In a preferred embodiment, the present invention may provide a nucleic acid detection system comprising the amplified nucleic acid preparations immobilized to a solid support.

25           In yet another preferred embodiment, the invention relates to a kit comprising a means for the single-phase amplification. Preferably, the means is a reaction vessel containing one or more reagents in concentrated form, where the reagent may be an enzyme or enzyme mix.

          Other objectives, features, and advantages of the present invention will  
30 become apparent from the following detailed description. The detailed description and the specific examples, while indicating preferred embodiments of the invention, are provided by way of illustration only. Accordingly, the present invention also

includes those various changes and modifications within the spirit and scope of the invention that may become apparent to those skilled in the art from this detailed description.

## 5 BRIEF DESCRIPTION OF THE FIGURES

Figure 1 depicts an overview of a preferred embodiment of the single-phase amplification methods of the present invention. PolyA<sup>+</sup> or total RNA is annealed with the single-stranded oligo dT-tailed promoter primer, T<sub>7</sub>T<sub>20</sub> (ggC Cag TgA ATT gTA ATA CgA CTC ACT ATA ggg Agg Cgg (T)<sub>20</sub>), creating a primer-template  
10 mixture. First strand cDNA synthesis is accomplished by combining the first strand cDNA reagent mix (Superscript II, buffer, DTT, and dNTPs) with the primer-template mixture and incubating at the appropriate time and temperature. Second strand cDNA synthesis is then performed by mixing the first strand cDNA reaction with second strand reagent mix, containing secondary cDNA mix (depc-H<sub>2</sub>O, Tris-HCl (pH7.0),  
15 MgCl<sub>2</sub>, (NH<sub>4</sub>)SO<sub>4</sub>, beta-NAD<sup>+</sup>, and dNTPs) and cDNA enzyme mix (Amplitaq DNA polymerase, *E. coli* ligase, *E. coli* RNase H, and *E. coli* DNA polymerase I)), followed by incubation at the appropriate times and temperatures. The resulting double-stranded (ds) cDNA contains a functional T7 RNA polymerase promoter, which is utilized for transcription. Finally, in vitro transcription is performed by  
20 combining the (ds) cDNA with IVT reagent (NTP, buffer, T7 RNA polymerase), yielding amplified, antisense RNA.

## DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENTS

The present invention can eliminate the need for multiple steps and varying  
25 reaction conditions and their associated problems. In preferred embodiments of the present invention, at least three otherwise separate enzymatic reactions can occur consecutively in one phase (*i.e.*, without organic extraction and precipitation), more preferably in the same reaction vessel. Preferably, cDNA synthesis according to the present invention may occur in a modified low salt buffer. In addition, the invention  
30 may involve a four-enzyme mix, which may include a thermal stable DNA polymerase and reverse transcriptase for the production of cDNA, terminal transferase, and RNA polymerase for RNA transcription. The enzyme activity may

be inactivated at the appropriate step with either heat or chemical treatment (*i.e.*, adjusting the salt concentration) or by the addition of an antibody specific to the enzyme. Furthermore, the single-phase amplification of nucleic acids can take place in the presence of less nucleotides than those which may be required in known methods.

Those skilled in the art will recognize that the products and methods embodied in the present invention may be applied to a variety of systems, including commercially available gene expression monitoring systems involving nucleic acid probe arrays, membrane blots, microwells, beads, and sample tubes, constructed with various materials using various methods known in the art. Accordingly, the present invention is not limited to any particular environment, and the following description of specific embodiments of the present invention are for illustrative purposes only.

In a preferred embodiment, the present invention can involve the amplification of nucleic acids, which may comprise synthesizing double-stranded DNA from a single-stranded DNA population, and producing multiple copies of RNA from the double-stranded DNA, where the synthesizing and producing occur in reaction vessels, preferably in the same reaction vessel.

The reaction vessel according to the present invention may include a membrane, filter, microscope slide, microwell, sample tube, array, or the like. *See* International Patent applications No. PCT/US95/07377 and PCT/US96/11147, which are expressly incorporated herein by reference. The reaction vessel may be made of various materials, including polystyrene, polycarbonate, plastics, glass, ceramic, stainless steel, or the like. The reaction vessel may preferably have a rigid or semi-rigid surface, and may preferably be conical (*e.g.*, sample tube) or substantially planar (*e.g.*, flat surface) with appropriate wells, raised regions, etched trenches, or the like. The reaction vessel may also include a gel or matrix in which nucleic acids may be embedded. *See* A. Mirzabekov *et al.*, *Anal. Biochem.* 259(1):34-41 (1998), which is expressly incorporated herein by reference.

Nucleic acids according to the present invention may include any polymer or oligomer of pyrimidine and purine bases, preferably cytosine, thymine, and uracil, and adenine and guanine, respectively. *See* Albert L. Lehninger, *Principles of Biochemistry*, at 793-800 (Worth Pub. 1982). Indeed, the present invention



contemplates any deoxyribonucleotide or ribonucleotide component, and any chemical variants thereof, such as methylated, hydroxymethylated or glucosylated forms of these bases, and the like. The polymers or oligomers may be heterogeneous or homogeneous in composition, and may be isolated from naturally-occurring sources or may be artificially or synthetically produced. *See* U.S. Patent No. 5,800,992, which is expressly incorporated herein by reference. In addition, the nucleic acids may be DNA or RNA, or a mixture thereof, and may exist permanently or transitionally in single-stranded or double-stranded form, including homoduplex, heteroduplex, and hybrid states.

The single-stranded or double-stranded DNA populations according to the present invention may refer to any mixture of two or more distinct species of single-stranded DNA or double-stranded DNA, which may include DNA representing genomic DNA, genes, gene fragments, oligonucleotides, polynucleotides, nucleic acids, PCR products, expressed sequence tags (ESTs), or nucleotide sequences corresponding to known or suspected single nucleotide polymorphisms (SNPs), having nucleotide sequences that may overlap in part or not at all when compared to one another. The species may be distinct based on any chemical or biological differences, including differences in base composition, order, length, or conformation. The single-stranded DNA population may be isolated or produced according to methods known in the art, and may include single-stranded cDNA produced from a mRNA template, single-stranded DNA isolated from double-stranded DNA, or single-stranded DNA synthesized as an oligonucleotide. The double-stranded DNA population may also be isolated according to methods known in the art, such as PCR, reverse transcription, and the like.

Where the nucleic acid sample contains RNA, the RNA may be total RNA, poly(A)<sup>+</sup> RNA, mRNA, rRNA, or tRNA, and may be isolated according to methods known in the art. *See, e.g.* T. Maniatis *et al.*, *Molecular Cloning: A Laboratory Manual*, 188-209 (Cold Spring Harbor Lab., Cold Spring Harbor, NY 1982, which is expressly incorporated herein by reference. The RNA may be heterogeneous, referring to any mixture of two or more distinct species of RNA. The species may be distinct based on any chemical or biological differences, including differences in base composition, length, or conformation. The RNA may contain full length mRNAs or

mRNA fragments (*i.e.*, less than full length) resulting from *in vivo*, *in situ*, or *in vitro* transcriptional events involving corresponding genes, gene fragments, or other DNA templates. In a preferred embodiment, the mRNA population of the present invention may contain single-stranded poly(A)+ RNA, which may be obtained from a RNA mixture (*e.g.*, a whole cell RNA preparation), for example, by affinity chromatography purification through an oligo-dT cellulose column.

Where the single-stranded DNA population of the present invention is cDNA produced from a mRNA population, it may be produced according to methods known in the art. *See, e.g.* Maniatis *et al.*, *supra*, at 213-46. In a preferred embodiment, a sample population of single-stranded poly(A)+ RNA may be used to produce corresponding cDNA in the presence of reverse transcriptase, oligo-dT primer(s) and dNTPs. Reverse transcriptase may be any enzyme that is capable of synthesizing a corresponding cDNA from an RNA template in the presence of the appropriate primers and nucleoside triphosphates. In a preferred embodiment, the reverse transcriptase may be from avian myeloblastosis virus (AMV), Moloney murine leukemia virus (MMuLV) or Rous Sarcoma Virus (RSV), for example, and may be thermal stable enzyme (*e.g.*, hTth DNA polymerase).

In a preferred embodiment of the present invention, the single-stranded cDNA produced using a mRNA population as template may be isolated from any resulting RNA:DNA heteroduplexes by heat or enzyme treatment (*e.g.*, RNase H). In a preferred embodiment, terminal transferase may be used to add poly(A) or poly(G) sequences to the 3'-termini of the single-stranded DNA. The double-stranded DNA of the present invention may be synthesized from the heterogeneous single-stranded DNA.

An oligonucleotide primer may be applied to the poly(A), poly(G), poly(C) or poly (T) tailed heterogeneous single-stranded DNA. The oligonucleotide primer preferably includes a poly(T) or poly(C) region complementary to the poly(A) or poly(G) tail attached to the single-stranded DNA. In addition, the oligonucleotide primer preferably includes a promoter consensus sequence capable of facilitating transcription by the RNA polymerase used, for example, the DNA-directed RNA polymerases derived from bacteriophage T7, T3 or SP6. The oligonucleotide primer may be synthesized, for example, using a PCR-MATE Model 391 DNA synthesizer

(Applied Biosystems) and purified by high-performance liquid chromatography before use. Second strand DNA synthesis may occur to yield the double-stranded DNA. *See, e.g., Examples, infra.*

5 In a preferred embodiment of the present invention, the ends of the double-stranded DNA may be blunted to prevent any concatenation of the double-stranded DNA. T4 DNA polymerase or *Escherichia coli* DNA polymerase I (Klenow fragment), for example, may be used preferably to produce blunt ends in the presence of the appropriate dNTPs.

10 In another preferred embodiment, multiple copies of the DNA may be obtained according to PCR methods known in the art in the presence of the appropriate primers. *See* R.K. Saiki *et al.*, *Science* 220:1350-1354 (1985), which is expressly incorporated herein by reference. In such circumstances, PCR cycles may preferably be limited to less than twenty to minimize amplification bias.

15 Multiple copies of RNA according to the present invention may be obtained by *in vitro* transcription from the DNA preferably using T7 RNA polymerase in the presence of the appropriate nucleoside triphosphates.

20 In a preferred embodiment of the present invention, the multiple copies of RNA may be labeled by the incorporation of biotinylated, fluorescently labeled or radiolabeled CTP or UTP during the RNA synthesis. *See* U.S. Patent No. 8,800,992 and International Patent Application PCT/US96/14839, which is expressly incorporated herein by reference. Alternatively, labeling of the multiple copies of RNA may occur following the RNA synthesis via the attachment of a detectable label in the presence of terminal transferase. In a preferred embodiment of the present invention, the detectable label may be radioactive, fluorometric, enzymatic, or  
25 colorimetric, or a substrate for detection (*e.g.*, biotin). Other detection methods, involving characteristics such as scattering, IR, polarization, mass, and charge changes, may also be within the scope of the present invention.

30 In a preferred embodiment, the amplified DNA or RNA of the present invention may be analyzed with a gene expression monitoring system. Several such systems are known. *See, e.g.*, U.S. Patent No. 5,677,195; Wodicka *et al.*, *Nature Biotechnology* 15:1359-1367 (1997); Lockhart *et al.*, *Nature Biotechnology* 14:1675-1680 (1996), which are expressly incorporated herein by reference. A gene

expression monitoring system according to the present invention may be a nucleic acid probe array such as the GeneChip® nucleic acid probe array (Affymetrix, Santa Clara, California). *See* U.S. Patents No. 5,744,305, 5,445,934, 5,800,992 and International Patent applications PCT/US95/07377, PCT/US96/14839, and

- 5 PCT/US96/14839, which are expressly incorporated herein by reference. A nucleic acid probe array preferably comprises nucleic acids bound to a substrate in known locations. In other embodiments, the system may include a solid support or substrate, such as a membrane, filter, microscope slide, microwell, sample tube, bead, bead array, or the like. The solid support may be made of various materials, including
- 10 paper, cellulose, nylon, polystyrene, polycarbonate, plastics, glass, ceramic, stainless steel, or the like. The solid support may preferably have a rigid or semi-rigid surface, and may preferably be spherical (*e.g.*, bead) or substantially planar (*e.g.*, flat surface) with appropriate wells, raised regions, etched trenches, or the like. The solid support may also include a gel or matrix in which nucleic acids may be embedded. *See* A.
- 15 Mirzabekov *et al.*, *Anal. Biochem.* 259(1):34-41 (1998).

The gene expression monitoring system, in a preferred embodiment, may comprise a nucleic acid probe array (including an oligonucleotide array, a cDNA array, a spotted array, and the like), membrane blot (such as used in hybridization analysis such as Northern, Southern, dot, and the like), or microwells, sample tubes,

20 beads or fibers (or any solid support comprising bound nucleic acids). *See* U.S. Patents No. 5,770,722, No. 5,744,305, No. 5,677,195 and No. 5,445,934, which are incorporated here in their entirety by reference. *See also* Examples, *infra*. The gene expression monitoring system may also comprise nucleic acid probes in solution.

The gene expression monitoring system according to the present invention

25 may be used to facilitate a comparative analysis of expression in different cells or tissues, different subpopulations of the same cells or tissues, different physiological states of the same cells or tissue, different developmental stages of the same cells or tissue, or different cell populations of the same tissue. *See* U.S. Patent No. 5,800,922. In a preferred embodiment, the proportional amplification methods of the present

30 invention can provide reproducible results (*i.e.*, within statistically significant margins of error or degrees of confidence) sufficient to facilitate the measurement of quantitative as well as qualitative differences in the tested samples. The proportional

amplification methods of the present invention may also facilitate the identification of single nucleotide polymorphisms (SNPs) (*i.e.*, point mutations that can serve, for example, as markers in the study of genetically inherited diseases) and other genotyping methods from limited sources. *See, e.g.*, Collins *et al.*, 282 *Science* 682 (1998), which is expressly incorporated herein by reference. The mapping of SNPs can occur by any of various methods known in the art, one such method being described in U.S. Patent No. 5,679,524, which is hereby incorporated by reference.

The RNA, single-stranded DNA, or double-stranded DNA population of the present invention may be obtained or derived from any tissue or cell source. Indeed, the nucleic acid sought to be amplified may be obtained from any biological or environmental source, including plant, viron, bacteria, fungi, or algae, from any sample, including body fluid or soil. In one embodiment, eukaryotic tissue is preferred, and in another, mammalian tissue is preferred, and in yet another, human tissue is preferred. The tissue or cell source may include a tissue biopsy sample, a cell sorted population, cell culture, or a single cell. In a preferred embodiment, the tissue source may include brain, liver, heart, kidney, lung, spleen, retina, bone, lymph node, endocrine gland, reproductive organ, blood, nerve, vascular tissue, and olfactory epithelium. In yet another preferred embodiment, the tissue or cell source may be embryonic or tumorigenic.

Tumorigenic tissue according to the present invention may include tissue associated with malignant and pre-neoplastic conditions, not limited to the following: acute lymphocytic leukemia, acute myelocytic leukemia, myeloblastic leukemia, promyelocytic leukemia, myelomonocytic leukemia, monocytic leukemia, erythroleukemia, chronic myelocytic (granulocytic) leukemia, chronic lymphocytic leukemia, polycythemia vera, lymphoma, Hodgkin's disease, non-Hodgkin's disease, multiple myeloma, Waldenstrom's macroglobulinemia, heavy chain disease, solid tumors, fibrosarcoma, myxosarcoma, liposarcoma, chondrosarcoma, osteogenic sarcoma, chordoma, angiosarcoma, endotheliosarcoma, lymphangiosarcoma, lymphangioendotheliosarcoma, synovioma, mesothelioma, Ewing's tumor, leiomyosarcoma, rhabdomyosarcoma, colon carcinoma, pancreatic cancer, breast cancer, ovarian cancer, prostate cancer, squamous cell carcinoma, basal cell carcinoma, adenocarcinoma, sweat gland carcinoma, sebaceous gland carcinoma,

papillary carcinoma, papillary adenocarcinomas, cystadenocarcinoma medullary carcinoma, bronchogenic carcinoma, renal cell carcinoma, hepatoma, bile duct carcinoma, choriocarcinoma, seminoma, embryonal carcinoma, Wilms' tumor, cervical cancer, testicular tumor, lung carcinoma, small cell lung carcinoma, bladder carcinoma, epithelial carcinoma, glioma, astrocytoma, medulloblastoma, craniopharyngioma, ependymoma, pinealoma, hemangioblastoma, acoustic neuroma, oligodendroglioma, menangioma, melanoma, neuroblastoma, and retinoblastoma. *See* Fishman *et al.*, *Medicine*, 2d Ed. (J.B. Lippincott Co., Philadelphia, PA 1985), which is expressly incorporated herein by reference.

10 In yet another preferred embodiment of the present invention, a nucleic acid detection system, the proportionally amplified DNA or RNA, or fragments thereof, may be immobilized directly or indirectly to a solid support or substrate by methods known in the art (*e.g.*, by chemical or photoreactive interaction, or a combination thereof). *See* U.S. Patent No. 5,800,992. The resulting immobilized RNA may be  
15 used as probes to detect nucleic acids in a sample population that can hybridize under desired stringency conditions. Such nucleic acids may include DNA contained in the clones and vectors of cDNA libraries.

The materials for use in the present invention are ideally suited for the preparation of a kit suitable for the single-phase proportional amplification of nucleic  
20 acids. Such a kit may comprise reaction vessels, each with one or more of the various reagents, preferably in concentrated form, utilized in the methods. The reagents may comprise, but are not limited to the following: low modified salt buffer, appropriate nucleotide triphosphates (*e.g.* dATP, dCTP, dGTP, dTTP; or rATP, rCTP, rGTP, and UTP) reverse transcriptase, RNase H, thermal stable DNA polymerase, RNA  
25 polymerase, and the appropriate primer complexes. In addition, the reaction vessels in the kit may comprise 0.2-1.0 ml tubes capable of fitting a standard PCR thermocycler, which may be available singly, in strips of 8, 12, 24, 48, or 96 well plates depending on the quantity of reactions desired. Hence, the single-phase amplification of nucleic acids may be automated, *e.g.*, performed in a PCR thermocycler. The PCR  
30 thermocyclers may include, but are not limited to the following: Perkin Elmer 9600, MJ Research PTC 200, Techne Gene E, Erichrom, and Whatman Biometra T1 Thermocycler.

Also, the automated machine of the present invention may include an integrated reaction device and a robotic delivery system. In such cases, part of all of the operation steps may automatically be done in an automated cartridge.

Without further elaboration, one skilled in the art with the preceding description can utilize the present invention to its fullest extent. The following examples are illustrative only, and not intended to limit the remainder of the disclosure in any way.

#### EXAMPLE ONE

##### Amplified Sample Preparation.

##### *cDNA synthesis*

Step 1: Primer-template annealing. The HPLC purified primer may be obtained from a -20°C storage stock, prepared in a 100 µM solution with TE (10 mM Tris-HCl (pH 8.0) and 1 mM EDTA (pH 8.0)) and diluted 1:1 with glycerol (for a final concentration of 50 µM in 50% glycerol and 50% TE). Where the desired nucleic acid sample is poly(A)+ RNA, a T<sub>7</sub>T<sub>20</sub> primer (ggc cag tga att gta ata cga ctc act ata ggg agg cgg (T)<sub>20</sub>) (Operon Technologies, Inc., Alameda, California) for example, may be used. In such case, an RNA sample (10 to 100 ng mRNA or 1-2 ug total RNA suspended in 2.5 µl or less) can be mixed with 0.5 µl primer to give a final volume of 3 µl. The mixture can be incubated at 70°C for 5-10 minutes, then cooled to 4°C.

Step 2: First strand cDNA synthesis. To the 3 µl primer-template mixture, may be added 2 µl of RT reagent, which may contain 0.5 µl Superscript II (200 U/µl, 100 U/reaction (Life Technologies, Inc., Gaithersburg, Maryland)) and 1.5 µl RT mix (100 µl 5X 1° cDNA buffer (shipped with Superscript II); 25 µl DTT (0.1 M); 25 µl dNTP (10 mM)). This 5 µl reaction mixture may then be incubated at 37°C for 60 minutes, then cooled to 4°C or on ice.

Step 3: Second strand cDNA synthesis. The 5 µl first strand cDNA reaction mixture may be mixed with 5 µl reagent mix (4.5 µl secondary cDNA mix and with 0.5 µl cDNA enzyme mix) at 4°C or on ice. Ensure that no more than 0.5 µl of enzyme is used. The reaction tube may be placed in a PCR cycler set up with the following program: extension for 2 hours at 16°C, then heat inactivation of enzyme

for 20 minutes at 75°C, and an ending temperature of 25°C. The extension step may be within the range of 10-37°C for many hours and the heat inactivation step may be within the range of 55-85°C for 1 minute to many hours.

The secondary cDNA mix may be prepared in 4.5 µl aliquots for 100 reactions by combining 245 µl RNase-free H<sub>2</sub>O, 40 µl of 1 M Tris HCl (pH 7.0), 35 µl of 0.1 M MgCl<sub>2</sub>, 100 µl of 0.1 M (NH<sub>4</sub>)<sub>2</sub>SO<sub>4</sub>, 15 µl of 10 mM beta-NAD<sup>+</sup>, and 15 µl of 10 mM dNTP.

The cDNA enzyme mix may be prepared in 0.5 µl aliquots for 100 reactions by combining in a screw-capped tube cooled to -20°C, 5 µl of Amplitaq DNA polymerase (5 U/µl) (PE Biosystems, Foster City, California), 5 µl of *E. coli* ligase (10 U/µl) (available from, for example, NEB, Inc., Beverly, Massachusetts), 6 µl of *E. coli* RNase H (2 U/µl) (available from, for example, Promega, Inc., Madison, Wisconsin), and 25 µl of *E. coli* DNA polymerase I (10 U/µl) (available from, for example, NEB, Inc., Beverly, Massachusetts). This mixture may be total volume adjusted to 50 µl with 9 µl of 50% cold glycerol, followed by a brief and gentle mix and quick spin before storage at -20°C.

#### *cRNA synthesis*

Step 4: *In vitro* transcription. To achieve maximal amplification, the 10 µl total volume of double-stranded cDNA may be combined with 10 µl of IVT reagent, which may contain 7 µl of NTP mix, 1 µl of 10X MEGAscript buffer (Ambion MEGAscript Kit, Ambion, Austin, Texas), and 2 µl of T<sub>7</sub> enzyme mix (Ambion MEGAscript Kit, Ambion, Austin, Texas). Alternatively, a dilution may be appropriate, for example, diluting 5 µl of double-stranded cDNA with 5 µl of H<sub>2</sub>O, and combining with 10 µl of IVT reagent.

The NTP mix may be prepared in 7 µl aliquots for 10 reactions, by combining 17.5 µl of 10 mM bCTP, 17.5 µl of 10 mM bUTP, 7.5 µl of 75 mM CTP, 7.5 µl of 75 mM UTP, 10 µl of 75 mM ATP, and 10 µl of 75 mM GTP.

The 20 µl total volume mixture may be incubated at 37°C for 4-6 hours, but can be 30 minutes to many hours. The resulting samples may be stored at -20°C, or analyzed.



Analysis may occur through the resolution of a 0.5  $\mu$ l or 1  $\mu$ l sample on a 1% agarose gel. Purification or quantification of the nucleic acid sample may occur by any one of the methods known in the art.

#### EXAMPLE TWO

##### 5 GeneChip® Analysis.

GeneChip® nucleic acid probe arrays are manufactured using technology that combines photolithographic methods and combinatorial chemistry. In a preferred embodiment, over 280,000 different oligonucleotide probes are synthesized in a 1.28 cm x 1.28 cm area on each array. Each probe type is located in a specific area on the probe array called a probe cell. Measuring approximately 24  $\mu$ m x 24  $\mu$ m, each probe  
10 cell contains more than  $10^7$  copies of a given oligonucleotide probe.

Probe arrays are manufactured in a series of cycles. A glass substrate is coated with linkers containing photolabile protecting groups. Then, a mask is applied that exposes selected portions of the probe array to ultraviolet light. Illumination  
15 removes the photolabile protecting groups enabling selective nucleotide phosphoramidite addition only at the previously exposed sites. Next, a different mask is applied and the cycle of illumination and chemical coupling is performed again. By repeating this cycle, a specific set of oligonucleotide probes is synthesized, with each probe type in a known physical location. The completed probe arrays are packaged  
20 into cartridges.

During the laboratory procedure, biotin-labeled RNA fragments referred to as the RNA target are hybridized to the probe array. The hybridized probe array is stained with streptavidin phycoerythrin conjugate and scanned by the Hewlett-Packard (HP) GeneArray™ Scanner at the excitation wavelength of 488 nm. The  
25 amount of emitted light at 570 nm and above is proportional to the amount of bound labeled target at each location on the probe array.

##### Step 1: Target Preparation

A total RNA population may be isolated from tissue or cells and reverse transcribed to produce cDNA. Then, *in vitro* transcription (IVT) produces biotin-  
30 labeled cRNA from the cDNA. The cRNA may be fragmented before hybridization.

## Step 2: Target Hybridization

After the biotin-labeled cRNA is fragmented, a hybridization cocktail is prepared, which includes labeled sample (0.05  $\mu\text{g}/\mu\text{l}$ ), probe array controls (1.5, 5, 25 and 100 pM respectively), herring sperm DNA (0.1 mg/ml), and BSA (0.5mg/ml). A cleanup procedure is performed on the hybridization cocktail after which 200  $\mu\text{l}$  is applied to the probe array through one of the septa in the array. It is then hybridized to the probes on the probe array during a 16-hour incubation at 45°C.

The hybridization protocol involves the following: (1) equilibrate probe array to room temperature immediately before use; (2) heat the sample(s) to 95°C for 5 minutes in a heat block; (3) meanwhile, wet the array by filling it through one of the septa with 1X Hybridization Buffer (1M NaCl, 0.1 M MES pH 6.7, 0.01% Triton X-100) using a micropipettor and appropriate tips; incubate the probe array at the hybridization temperature for 10 minutes with rotation; (5) after incubation at 95°C (step #2 above), transfer the samples to a 45°C heat block for 5 minutes; (5) spin samples at maximum speed in a microcentrifuge for 5 minutes to remove any insoluble material from the hybridization mixture; (6) remove the buffer solution from the probe array cartridge and fill with 200  $\mu\text{l}$  of the clarified hybridization cocktail avoiding any insoluble matter in the 20  $\mu\text{l}$  at the bottom of the tube; (7) place probe array in rotisserie box in 45°C oven; load probe arrays in a balanced configuration around rotisserie axis; rotate at 60 rpm; and (8) hybridize for 16 to 40 hours.

## Step 3: Probe Array Washing, Staining, and Fluidics Station Setup

Immediately following the hybridization, the hybridized probe array undergoes manual washing and staining, then washing on the fluidics station. The protocol involves the following: (1) remove the hybridization cocktail from the probe array and set it aside in a microcentrifuge tube; store on ice during the procedure or at -20°C for long-term storage; (2) rinse the probe array by pipetting 200  $\mu\text{l}$  1X MES buffer pH 6.7 through one of the probe array septa; (3) fill the probe array septa with 200  $\mu\text{l}$  6X SSPE-T (300 ml of 20X SSPE and 500  $\mu\text{l}$  of 10% Triton X 100 to 700 ml of water, final pH 7.6) and wash with 6X SSPE-T on the fluidics station with wash A cycle (10 cycles, drain and fill twice each cycle); (4) remove the 6X SSPE-T and rinse the probe array with 0.1X MES buffer pH 6.7 (0.1 M MES, 0.1 M NaCl and 0.01%

Triton); (5) fill the probe array with 200 µl 0.1X MES and incubate at 45°C on the rotisserie at 60 rpm for 30 minutes; and (6) remove the 0.1X MES, rinse the probe array with 1X MES in the probe array while preparing the stain.

- 5 Staining the probe array involves preparing Streptavidin Phycoerythrin (SAPE) stain solution. Stain should be stored in the dark and foil wrapped or kept in an amber tube at 4°C. Remove stain from refrigerator and tap the tube to mix well before preparing stain solution. The concentrated stain or diluted SAPE stain solution should not be frozen. The SAPE stain should be prepared immediately before use.

- 10 For each probe array to be stained, combine the following components to a total volume of 200 µl (1:100 dilution of SAPE, final concentration of 10 µg/ml): 188 µl 1X MES; 10 µl of 50 mg/ml acetylated BSA (final concentration of 2.5 mg/ml); and 2 µl of 1 mg/ml streptavidin phycoerythrin (SAPE).

Remove the 1X MES and apply the stain solution to the probe array. Incubate for 15 minutes at 60 rpm at room temperature or 40°C.

- 15 Remove the stain and fill the probe array with 6X SSPE-T. Wash the probe array with 6X SSPE-T on the fluidics station with wash A cycle.

- The experiment parameters are preferably defined using commercially available GeneChip® software (Affymetrix, Santa Clara, California) on a PC-compatible workstation with a Windows NT® operating system. The probe array type, sample description, and comments are entered in the software and saved with a unique experiment name.

- The user protocol involves the following: (1) launch the software from the workstation and choose Experiment Info from the Run menu; alternatively, click the New Experiment icon on the GeneChip® software tool bar; the Experiment Information dialog box will appear allowing the experiment name to be defined along with several other parameters such as probe array type, sample description, and comments; (2) type in the experiment name; click on the box to the right of Probe Array type and select the probe array type from the drop-down list; experiment name and probe array type are required; complete as much of the other information as desired; the protocol information at the bottom of the dialog box will be imported to the experiment information dialog box after the hybridization and scan have been

completed; (3) save the experiment by choosing Save; the name of the experiment will be used by the software to access the probe array type and data for the sample while it is being processed; data files generated for the sample will be automatically labeled to correspond to the experiment name; the Protocol section of the dialog box will be filled in by the software; and (4) close the Experiment Information dialog box.

The GeneChip® Fluidics Station 400 is preferably used to wash the probe arrays. It is operated using the GeneChip® software as follows: (1) choose Fluidics from the Run menu; alternatively, click the Start Protocol icon on the GeneChip® software tool bar; the Fluidics Station dialog box will appear with a drop-down list for the experiment name; a second list is accessed for the Protocol for each of the four fluidics station modules; (2) prime the fluidics station, by clicking Protocol in the Fluidics Station dialog box; choose Prime for the respective modules in the Protocol drop-down list; change the intake buffer reservoir A and B to 6X SSPE-T; click Run for each module to begin priming; priming should be done whenever the fluidics station is first started up, when wash solutions are changed, after washing if a shutdown has been performed on any module, and if the LCD window instructs the user to prime; priming ensures that the wash lines are filled with the appropriate buffer and the fluidics station is ready for washing; a prime takes approximately 3 to 5 minutes to complete; the fluidics station LCD window and the Fluidics Station dialog box will display the status of the prime and give instructions as it progresses; follow the instructions on the LCD window and dialog box; when priming is complete, the LCD window and dialog box will indicate that the fluidics station is ready to run a wash; (3) wash the probe array on the fluidics station, by customizing the HYBWASH protocol to create a wash of 10 cycles with 2 mixes per cycle with 6X SSPE-T at room temperature; in the Fluidics Station dialog box on the workstation, select the correct experiment name in the drop-down Experiment list; the probe array type will appear automatically; in the Protocol drop-down list, select the modified HYBWASH protocol created in step 1 to control the wash of the probe array; if a customized protocol is run, check the parameters of each of the protocols chosen to be sure they are appropriate for your experiment; this can be done in the Fluidics Protocol dialog box found by choosing Edit Protocol under the Tools menu; choose Run in the Fluidics Station dialog box to begin the wash; follow the instructions on

- the LCD window on the fluidics station; open the probe array holder by pressing down on the probe array lever to the Eject position; place the appropriate probe array into the probe array holder of the selected module and gently push up on the lever to engage it; the latch should be secure when the probe array holder is fully closed; a
- 5 light click should be heard; engage the probe array holder lever by firmly pushing up on it to the Engage position; the Fluidics Station dialog box and the LCD window will display the status of the wash as it progresses; when the wash is complete, the LCD window will display EJECT CARTRIDGE; eject the probe array by pushing down firmly on the probe array lever; and (4) perform the cleanout procedure, by returning
- 10 the probe array to the probe array holder; latch the probe array holder by gently pushing it up until a light click is heard; engage by firmly pushing up on the probe array lever to the Engage position; the fluidics station will drain the probe array and then fill it with a fresh volume of the last wash buffer used; when it is finished, if the LCD window displays EJECT CARTRIDGE again, remove the probe array and
- 15 inspect it again for bubbles; if no bubbles are present, it is ready to scan; after ejecting the probe array from the probe array holder, the LCD window will display ENGAGE WASHBLOCK; latch the probe array holder by gently pushing it up and in until a light click is heard; engage the washblock by firmly pushing up on the probe array lever to the Engage position; the fluidics station will automatically perform a
- 20 Cleanout procedure; the LCD window will indicate the progress of the Cleanout procedure; when the Cleanout procedure is complete, the LCD window should display Washing done, READY; if no other washes are to be performed, place wash lines into a bottle filled with deionized water; choose Shutdown for all modules from the drop-down Protocol list in the Fluidics Station dialog box; click the Run button for all
- 25 modules; after Shutdown protocol is complete, flip the ON/OFF switch of the fluidics station to the OFF position; and scan the probe array.

#### Step 4: Probe Array Scan

- Once the probe array has been hybridized, stained, and washed, it is scanned. Each workstation running the software can control one scanner. Each scan takes
- 30 approximately 5 minutes, and two scans are recommended.

The scanner acquires an image of each of the hybridized 24  $\mu\text{m}$  x 24  $\mu\text{m}$  probe cells. Each complete probe array image is stored in a separate data file that

corresponds to its experiment name and is saved with a data image file (.dat) extension.

The scanner is also controlled by the GeneChip® software. The probe array is scanned after the wash protocols are complete. The probe array scan proceeds as follows: (1) choose Scanner from the Run menu; alternatively, click the Start Scan icon in the GeneChip® software tool bar; the Scanner dialog box will appear with a drop-down list of experiments that have not been run; a scrollable window will also be displayed showing previous scans; choose the experiment name that corresponds to the probe array to be scanned; a previously run experiment can also be chosen from the Previous Experiments list by double-clicking on the name desired; (2) check for the correct pixel value and wavelength of the laser beam; for a 24 µm x 24 µm probe array with a phycoerythrin stain: Pixel value = 3 µm, Wavelength = 570 nm; (3) once the experiment has been selected, click the Start button; a dialog box will prompt the user to load a sample into the scanner; and (4) load the Probe Array into the HP GeneArray™ Scanner; open the sample door on the scanner and insert the probe array into the holder; do not force the probe array into the holder; close the sample door of the scanner; start the Scan, by clicking OK in the Start Scanner dialog box; the scanner will begin scanning the probe array and acquiring data; when Scan in Progress is chosen from the View menu, the probe array image will appear on the screen as the scan progresses.

#### Step 5: Data Analysis and Interpretation

Data is analyzed using GeneChip® software. In the Image window, a grid is automatically placed over the image of the scanned probe array to demarcate the probe cells. After grid alignment (the user may adjust the alignment if necessary), the mean intensity at each probe cell is calculated by the software. The intensity patterns are analyzed.

After scanning the probe array, the resulting image data created is stored on the hard drive of the GeneChip® workstation as a .dat file with the name of the scanned experiment. In the first step of the analysis, a grid is automatically placed over the .dat file so that it demarcates each probe cell. One of the probe array library files, the .cif file, indicates to the software what size of grid should be used. Confirm

the alignment of the grid by zooming in on each of the four corners and the center of the image.

If the grid is not aligned correctly, adjust its alignment by placing the cursor on an outside edge or corner of the grid. The cursor image will change to a small double-headed arrow. The grid can then be moved using the arrow keys or by clicking and dragging its borders with the mouse.

Sample analysis occurs as follows: (1) choose Defaults from the Tools menu to access the Probe Array Call Settings tab dialog box; in the Defaults dialog box, click on the Probe Array Call Settings tab to display probe array calling algorithm choices; (2) highlight GeneChip® Expression and click the Modify button or double click the algorithm name; (3) in the Probe Array Call Settings dialog box, select the probe array type in the drop down list; for that probe array make sure the Use As Current Algorithm check box is selected; (4) click the OK button to apply your choices for the selected probe array type; (5) in the Defaults dialog box, click the OK button to apply your choices regarding parameters set by all of the tab dialog boxes in the window; (6) after confirming that the above parameters are correct, select the appropriate image to be analyzed; and (7) select Analysis from the Run menu or click the Run Analysis icon on the GeneChip® software tool bar; the software calculates the average intensity of each probe cell using the intensities of the pixels contained in the cell; pixels on the edges of each cell are not included, which prevents neighboring cell data from affecting a cells calculated average intensity; the calculated average intensity is assigned an X/Y-coordinate position, which corresponds to the cell's position on the array; this data is stored as a .cel file using the same name as the .exp and .dat files; the .cel file is an intermediate data file; the software then applies the selected probe array algorithm to determine expression levels for each gene; this is done with reference to the information contained in the .cdf file, the second library file for the probe array; the resulting analysis is automatically displayed as a .chp file in the Expression Analysis window of GeneChip® software; the .chp file has the same name as the .exp, .dat, and .cel files.

The specific embodiments described above do not limit the scope of the present invention in any way as they are single illustrations of individual aspects of the invention. Functionally equivalent methods and components are within the scope

of the invention. The scope of the appended claims thus includes modifications that will become apparent to those skilled in the art from the foregoing description.



We claim:

1. A method for the amplification of nucleic acid, said method comprising:  
synthesizing double-stranded DNA from a single-stranded  
DNA population, and  
producing multiple copies of RNA from said double-stranded  
DNA,  
wherein said amplification occurs in a single phase.
2. The method of claim 1, wherein said amplification is proportional.
3. The method of claim 1, wherein said nucleic acid is selected from the  
group consisting of genomic DNA, cDNA, total RNA, poly(A)<sup>+</sup> RNA, and  
oligonucleotides.
4. The method of claim 3, wherein said poly(A)<sup>+</sup> RNA is mRNA.
5. The method of claim 1, further comprising:  
contacting said multiple copies of RNA with a solid support  
comprising nucleic acid probes.
6. The method of claim 5, further comprising:  
detecting the presence or absence of hybridization of said multiple  
copies of RNA to said nucleic acid probes on said solid support.
7. The method of claim 5, wherein said solid support comprising nucleic acid  
probes is selected from the group consisting of a nucleic acid probe array, a  
membrane blot, a microwell, a bead, and a sample tube.
8. The method of claim 1, wherein said nucleic acid is isolated from an  
eukaryotic cell or tissue.

9. The method of claim 8, wherein said eukaryotic cell or tissue is mammalian.

10. The method of claim 9, wherein said mammalian cell or tissue is human.

11. The method of claim 1, wherein said nucleic acid is isolated from a source selected from the group consisting of dissected tissue, microdissected tissue, a tissue subregion, a tissue biopsy sample, a cell sorted population, a cell culture, and a single cell.

12. The method of claim 1, wherein said nucleic acid is isolated from a cell or tissue source selected from the group consisting of brain, liver, heart, kidney, lung, spleen, retina, bone, lymph node, endocrine gland, reproductive organ, blood, nerve, vascular tissue, and olfactory epithelium.

13. The method of claim 1, wherein said nucleic acid is isolated from a cell or tissue source selected from the group consisting of embryonic and tumorigenic.

14. An amplified nucleic acid preparation comprising RNA obtained by the method of claim 1.

15. An amplified nucleic acid preparation comprising RNA obtained by the method of claim 2.

16. A gene expression monitoring system comprising a solid support, which comprises nucleic acid probes and the amplified nucleic acid preparation of claim 14.

17. A gene expression monitoring system comprising a solid support, which comprises nucleic acid probes and the amplified nucleic acid preparation of claim 15.

18. A nucleic acid detection system comprising the amplified nucleic acid preparation of claim 14 immobilized to a solid support.

19. A nucleic acid detection system comprising the amplified nucleic acid preparation of claim 15 immobilized to a solid support.

20. The method of claim 1, wherein said synthesizing and producing comprise the use of an automated machine.

21. The method of claim 20, wherein said automated machine is selected from the group consisting of a PCR thermocycler, an integrated reaction device, and a robotic delivery system.

22. A kit for the amplification of nucleic acids, wherein said kit comprises a means for the single-phase amplification according to claim 1.

23. The kit of claim 22, wherein said means is a reaction vessel containing one or more reagents in concentrated form.

24. The kit of claim 23, wherein said reagent is an enzyme.

## ABSTRACT

The amplification of nucleic acids in a single phase can reliably provide products at relatively low cost and labor. A single-phase amplification can also increase the amount of nucleic acids while preserving the relative abundance of the individual nucleic acid species, or portions thereof. A single-phase amplified nucleic acid preparation may be analyzed in a gene expression monitoring system, preferably involving a nucleic acid probe array.

09373984 081699

# Single phase anti-sense RNA amplification

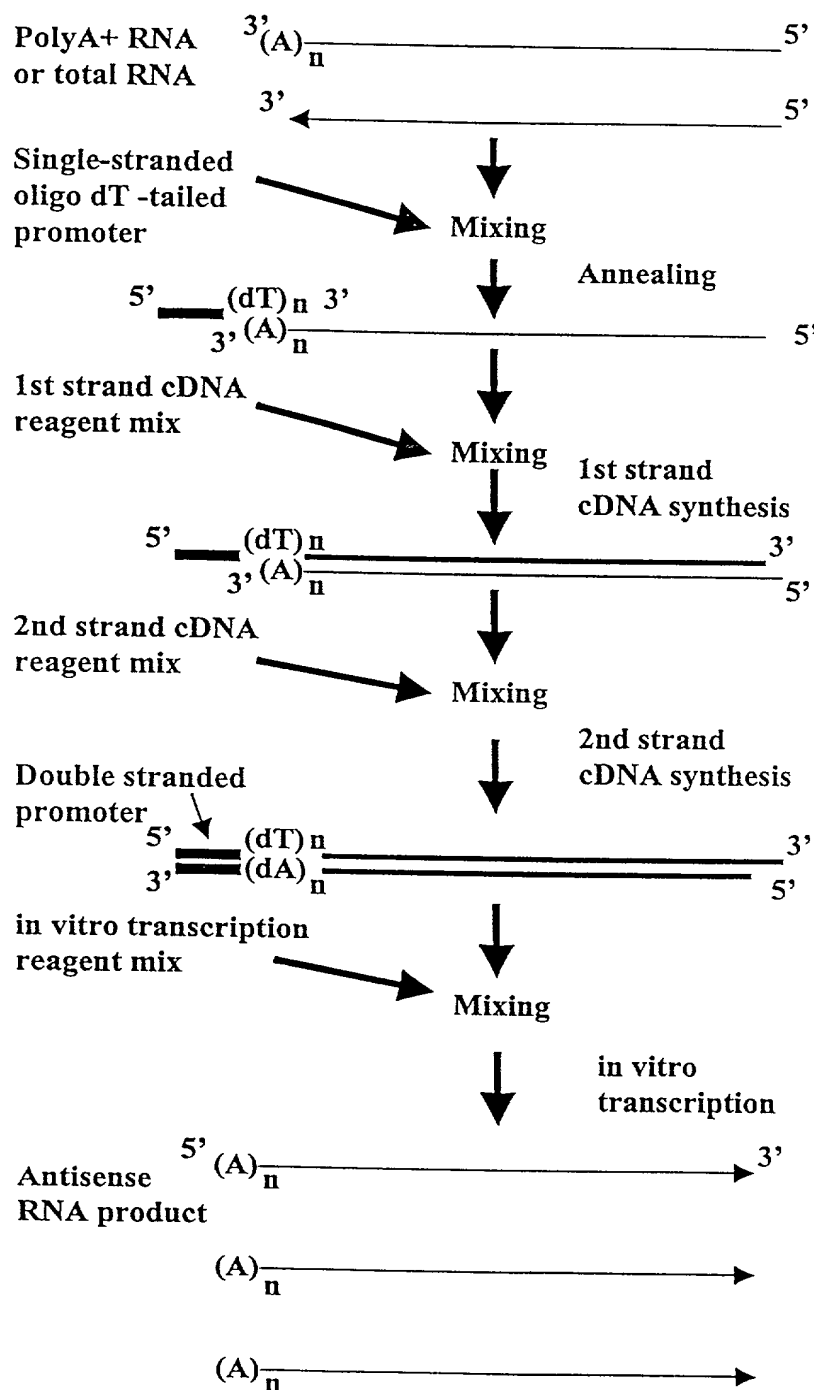


Figure 1

**DECLARATION AND POWER OF ATTORNEY**

As a below named inventor, I hereby declare that:

My residence, post office address and citizenship are as stated below next to my name.

I believe I am the original, first and sole inventor (if only one name is listed below) or an original, first and joint inventor (if plural names are listed below) of the subject matter which is claimed and for which a patent is sought on the invention entitled:

**SINGLE-PHASE AMPLIFICATION OF NUCLEIC ACIDS**

the specification of which is attached hereto unless the following box is checked:

☐ was filed on \_\_\_\_\_ as United States Application Number or PCT International Application Number \_\_\_\_\_ and was amended on \_\_\_\_\_ (if applicable).

I hereby state that I have reviewed and understand the contents of the above-identified specification, including the claims, as amended by any amendment referred to above.

I acknowledge the duty to disclose information which is known by me to be material to patentability as defined in Title 37, Code of Federal Regulations § 1.56.

I hereby claim foreign priority benefits under Title 35, United States Code, § 119(a)-(d) or § 365(b) of any foreign application(s) for patent or inventor's certificate, or § 365(a) of any PCT International application which designated at least one country other than the United States, listed below and have also identified below any foreign application for patent or inventor's certificate, or PCT International application having a filing date before that of the application on which priority is claimed:

**PRIOR FOREIGN APPLICATION(S)**

NUMBER	COUNTRY	DAY/MONTH/YEAR FILED	PRIORITY CLAIMED

I hereby claim the benefit under Title 35, United States Code § 119(e) of any United States provisional application(s) listed below.

APPLICATION NO.	FILING DATE
60/146,142	30 July 1999

0037394-03169  
669780-49627660


I hereby claim the benefit under Title 35, United States Code, § 120 of any United States application(s), or § 365(c) of any PCT International application designating the United States, listed below and, insofar as the subject matter of each of the claims of this application is not disclosed in the prior United States or PCT International application in the manner provided by the first paragraph of Title 35, United States Code, § 112, I acknowledge the duty to disclose information which is known by me to be material to patentability as defined in Title 37, Code of Federal Regulations § 1.56 which became available between the filing date of the prior application and the national or PCT International filing date of this application:

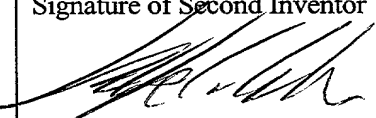
APPLICATION SERIAL NO.	FILING DATE	STATUS: PATENTED, PENDING, ABANDONED

I hereby appoint as my attorneys, with full powers of substitution and revocation, to prosecute this application and transact all business in the Patent and Trademark Office connected therewith: Don J. Pelto, Reg. No. 33,754, Jeff E. Schwartz, Reg. No. 39,019, and Lawrence M. Sung, Reg. No. 38,330.

Address all correspondence to Don J. Pelto, 1801 K Street, NW, Suite 400K, Washington, D.C. 20006-1301. Address telephone communications to Don J. Pelto at (202) 775-7100.

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under Section 1001 of Title 18 of the United States Code and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

Full Name of First or Sole Inventor Xing SU	Signature of First or Sole Inventor 	Date 8/12/99
Residence Address (City, State) Cupertino, CA		Country of Citizenship P.R. China
Post Office Address 891 E. Estates Drive, Cupertino, CA 95014		

Full Name of Second Inventor Rolfe C. ANDERSON	Signature of Second Inventor 	Date 8/11/99
Residence Address (City, State) Saratoga, CA		Country of Citizenship U.S.A.
Post Office Address 13690 Rossmere Court, Saratoga, CA 95020		

659130-4866260